



## **ATLANTIS**

**MB/PC250 - MB/PC260**

**INSTRUCTION MANUAL**

## WARNINGS

### Safety during installation

- Before installing and using the treatment chair, carefully read the instruction manual and ensure that all directions and instructions contained are fully understood. It includes important safety information for installation, operation and maintenance!
- When unpacked, check the contents of the delivery to ensure that all parts have been included and that nothing has been damaged during transit. In the event of any doubt, please contact our customer service department.
- The Atlantis treatment chair must be positioned on a flat surface, levelled and able to bear the weight when the chair is loaded, referring to the maximum allowable load.
- **Caution! Never move the Atlantis pedicure chair as shown in figure 3 but always as shown in figure 2.**

### Safety during use

- If the treatment chair is used for any purpose other than that for which it was designed and as described in this instruction manual, such use will invalidate the manufacturer's warranty and guarantee!
- This unit may only be used by personnel that complete the treatment from here on 'operator'.
- Children should not be left unattended in proximity to the equipment.

### Safety during maintenance

- Cleaning and maintenance must be made only by those qualified for such a job.

### Transport/ Storage/ Operating Requirements

Transport only in the manufacturer's original packaging.

Transport and storage temperature: -20 °C to +50 °C

Operating temperature: +10 °C to +40 °C

Relative air humidity: max. 90 %, non-condensing

To avoid the formation of condensation, abrupt changes of temperature should be avoided.

### FIELDS OF USE


The Atlantis treatment chair is a very comfortable pedicure chair for foot care and beauty treatments. Atlantis offers both practicality and style, with a large, fully upholstered chair that rotates 90° for ease of use for both the client and the manicure / pedicure specialist. The stylish chair is supported by a rectangular structure covered in stainless steel; also, on top of this elegant base is a beautifully designed basin in high quality, resistant, white techno polymer, complete with tap, shower head and floor-connected drainage system. The Atlantis also offers a padded, height-adjustable footrest that can also be moved closer or further away from the chair itself, to accommodate for clients of different heights, and - ONLY FOR VERSION MB/PC260 - a shiatsu massage system.

## COMPONENTS

### Legend

1. Stainless steel base
2. Plastic basin
3. Tap, shower head and floor-connected drainage system
4. Padded movable footrest
5. Rotating chair (90°)
6. Button panel for massage control (ONLY FOR VERSION MB/PC260)

Fig. 1

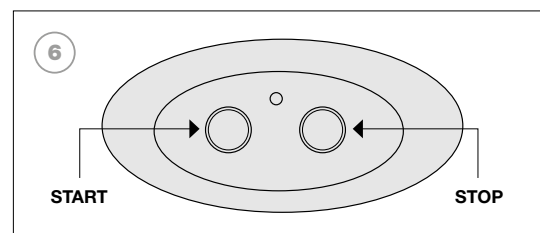


## INSTALLATION

**Caution! Never move the Atlantis pedicure chair as shown in figure 3 but always as shown in figure 2.**

Fig. 2



Fig. 3

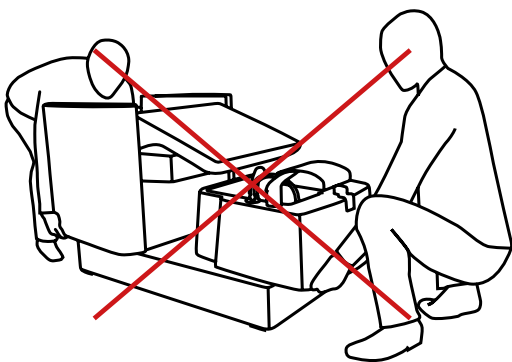


To install the pedicure chair:

1. Remove the chair's magnetised rear panel and mount the backrest
2. Remove the front / back steel panels on the base
3. Place the semi-rotation stopper on the desired side
4. Adjust the feet, if necessary
5. Trace the outline and fix the unit with the four adjustable angles
6. Connect the drain and water supply, see fig. 4
7. Check the plumbing
8. Restore the panelling

Fig. 4

- A= COLD WATER 3/8"
- B= DRAIN 0/40
- C= HOT WATER 3/8"
- D= ELECTRICAL SOCKET (ONLY FOR VERSION MB/PC260)

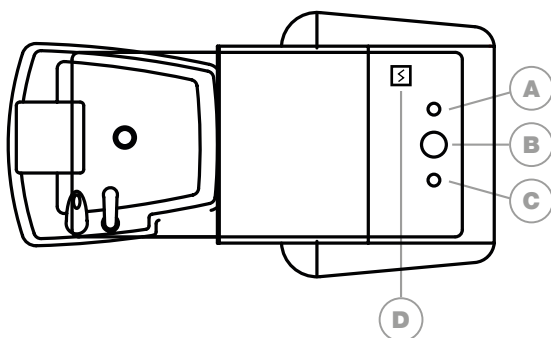

**USE**

1. Rotate the chair to make the customer sit
2. Regulate the height and distance of the footrest cushion
3. (ONLY FOR VERSION MB/PC260) If required, activate the massage as illustrated in the paragraph "Massage system Instructions"
4. Proceed to treatment

**MAINTENANCE**

**Cleaning and Maintenance**


Clean the upholstery of the treatment chair with warm water and a cloth.  
 Stains should be removed immediately to avoid them penetrating the material.  
 Clean the frame of the treatment chair with a lightly moistened cloth.  
 Do not rub the disinfectant in the upholstery.  
 Do not use any solvents, chlorides, polish, synthetic cleaning agents, wax polish, aerosol sprays or aggressive cleaning agents which could cause damage to the upholstery and the treatment chair.  
 Do not use any substances which contain oil or grease to clean the upholstery.  
 The upholstery is not resistant to chemical cleaning and dry-cleaning methods.

**Repairs and Servicing**

Periodically inspect and clean the plumbing.  
 The treatment chair does not require any special servicing when used correctly and in compliance with the manufacturer's directions and instructions.  
 (ONLY FOR VERSION MB/PC260) Every six months inspect the PVC cloth where the massage wheels flow and, depending on the frequency of use, lubricate it with synthetic rubber such as TOPAS NB 52, as shown in fig. 5.

Fig. 5

REMOVE THE CASING AND LUBRICATE THE PVC BACK CLOTH OF THE MESSAGE

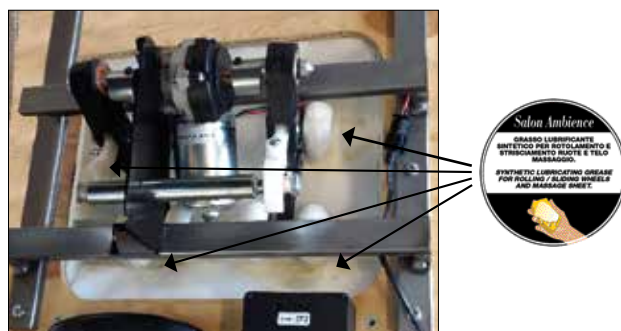

However, in the event of any problem arising, please contact our customer service department.  
 Only persons authorised by the manufacturer may perform repairs. Any alterations or repairs performed by non-authorised persons will invalidate any guarantee, warranty or responsibility of the manufacturer. In addition, the safety of the treatment chair could be adversely affected.

**Disposal**

Dispose of the packaging, treatment chair, upholstery and accessories in an environmentally friendly manner and in compliance with local regulations.

## ONLY FOR VERSION MB/PC260

### EQUIPMENT FEATURES



<b>Manufacturer</b>	Salon Habitat (Italia) s.r.l. via Enzo Lazzaretti 12 42122 Reggio Emilia www.medicalandbeauty.com info@medicalandbeauty.com
<b>Model</b>	Atlantis MB/PC260
<b>Function</b>	Treatment chair for foot care and beauty treatments equipped with shiatsu massage system
<b>IEC Protection Class</b>	II
<b>Year of Manufacture</b>	2020
<b>Power Supply</b>	100/220-240 V ~ 50-60 Hz
<b>Nominal Power</b>	110 W
<b>Maximum allowable load for the chair</b>	130 Kg

### POWER SUPPLY

The power supply of the Atlantis treatment chair has been created for an electric cable directly connected to the control unit, equipped with a plug for the socket.

Cable and plug are supplied with the treatment chair.

The power supply must be sized for the following characteristics – 100/220-240 V ~ 50-60 Hz, monophased. Furthermore, the electrical system must comply with the laws and regulations in force at the place of use, and must be equipped with a differential switch that ensures people's protection against indirect electrical contact. The power supply line of the Atlantis treatment chair must be protected against overload and short circuit. Plug the machine into an easily accessible socket, checking that the voltage (V) of the network corresponds to the plate on the machine.

- Avoid use of extension cords, adapters and multiple sockets; when possible insert the plug directly in the socket.
  - Keep the power cord away from areas of passage in order to avoid crushing, interference or obstacles to movement.
  - Keep the power cord away from heat sources or hot surfaces.
- To avoid dangerous electrical shocks, use a socket with correctly covered wires.

We ask that you request support from your trusted installation company.

## ELECTRICAL COMPONENTS



1. Power supply  
7333



2. Control-box  
3885




3. Massage system  
3888



4. Massage control panel  
3886

### MESSAGE SYSTEM INSTRUCTIONS

Fig. 6

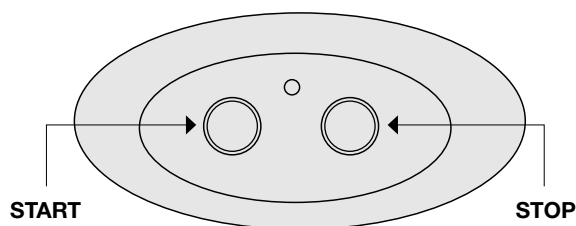

If requested, use the appropriate push-button panel to start and stop the massage movement, as shown in fig. 6.

Each button has to be held down to for 1 second IN ORDER TO AVOID UNWANTED ACTIONS.

The massage cycle lasts 3 minutes. To stop the massage, use the stop-button, as shown in fig. 6. AFTER 3 MINUTES IT WILL STOP INDEPENDENTLY.

### IF

- Strange noises such as buzzing or creaking
  - Strange smells such as burning
- come from the treatment chair

### OR

- Parts of the treatment chair get too warm during the normal functioning of equipment, immediately stop any activity, remove the plug from the socket and inform Salon Habitat Italia Srl Customer Service.



## DICHIARAZIONE DI CONFORMITA'

Noi: Salon Habitat (Italia) s.r.l.  
Via Lazzaretti 12  
42122 Reggio Emilia

Attestiamo sotto la nostra completa responsabilità che il prodotto:

**Poltrona podologica con massaggio Atlantis MB/P260**

cui la dichiarazione si riferisce risulta in conformità con quanto previsto dalle seguenti direttive comunitarie, comprese le ultime modifiche e con la relativa legislazione di recepimento:

- **Direttiva 2006/42/CE** \_  
concernente il ravvicinamento delle legislazioni degli stati membri relative alle macchine
- **Direttiva 2006/95/CE** \_  
concernente il ravvicinamento delle legislazioni degli stati membri relative al materiale elettrico destinato ad essere adoperato entro taluni limiti di tensione
- **Direttiva 2004/108/CE** \_  
concernente il ravvicinamento delle legislazioni degli stati membri relative alla compatibilità elettromagnetica e che abroga la direttiva 89/366/CEE

Reggio Emilia, 01/01/2022

Il legale rappresentante

Sig. Michael Frott  




**Salon Habitat (Italia) s.r.l. via Lazzaretti 12 42122 Reggio Emilia**  
**Sales Italia 0522 276740 International Sales 0522 276755 Fax 0522 276767**  
**info@medicalandbeauty.com • www.medicalandbeauty.com**